Clinical Case

CM Alvim NeoPoros

Initial Situation

Dr. Pedro Gomes

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upper

lower

Surgical products	Prosthetic products
CM Alvim NeoPoros	CM Universal abutment Exact
4.3 X 16mm	3.3 X 4 X 3.5mm
Xenograft	

- Profile: C. S.
Female 32 years old
Porto, Portugal
- Clinical Situation:
Internal resorption of tooth 11
- Restorative Solution:
Extraction, followed by immediate implant placement, immediate loading with a provisional crown on a universal abutment and bone substitute material. The case used “one-abutment-one-time” concept and finalized with a metal-ceramic crown.

Summary	Immediate loading protocol	
Surgical Description	Maxilla	Anterior
	1 Implants	Conventional
Restorative Solution	Conventional	Metal - ceramic

Clinical Case


Pre-operative X-ray

Tooth extraction

Implant and abutment in place

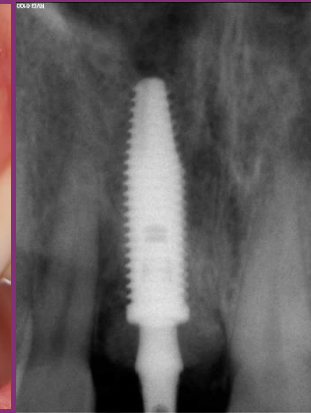
Post-operative X-ray

Xenograft material (Xenograft)

Immediate provisional restoration in place

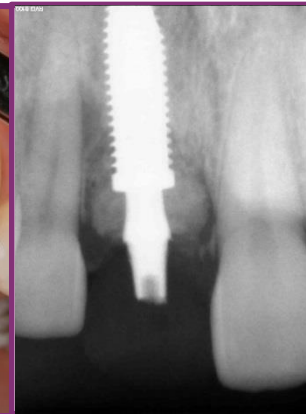Post provisional restoration X-ray

Provisional restoration 6-month after surgery

Clinical Case

*Provisional restoration 6-month after surgery
occlusal view*

*Stable soft tissue contour
6 months after surgery*


*Great emergence profile
6-month after surgery*

Final restoration 11 months after surgery

Final restoration

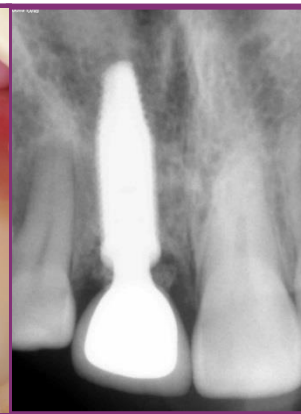*Final restoration
X-ray*

Extra oral frontal view

Clinical Case

4-year follow-up X-ray

4-year follow up - intraoral frontal view