

## INSTALLATION GUIDE FOR OPFLEX403 AND OPFLEX503 LIBRARY


### 3Shape®

#### Installation:

- Unzip the file: "Anthogyr-Labside-OPFLEX403\_503-SA-DT-PA-3Shape2018.zip" for 2018 3Shape version.  
"Anthogyr-Labside-OPFLEX403\_503-SA-DT-PA-3Shape2019.zip" for 2019 3Shape version.  
"Anthogyr-Labside-OPFLEX403\_503-SA-DT-PA-3Shape2020.zip" for 2020 or higher 3Shape version.

- Open 3Shape Dental System Control Panel.

- Go to import/export:

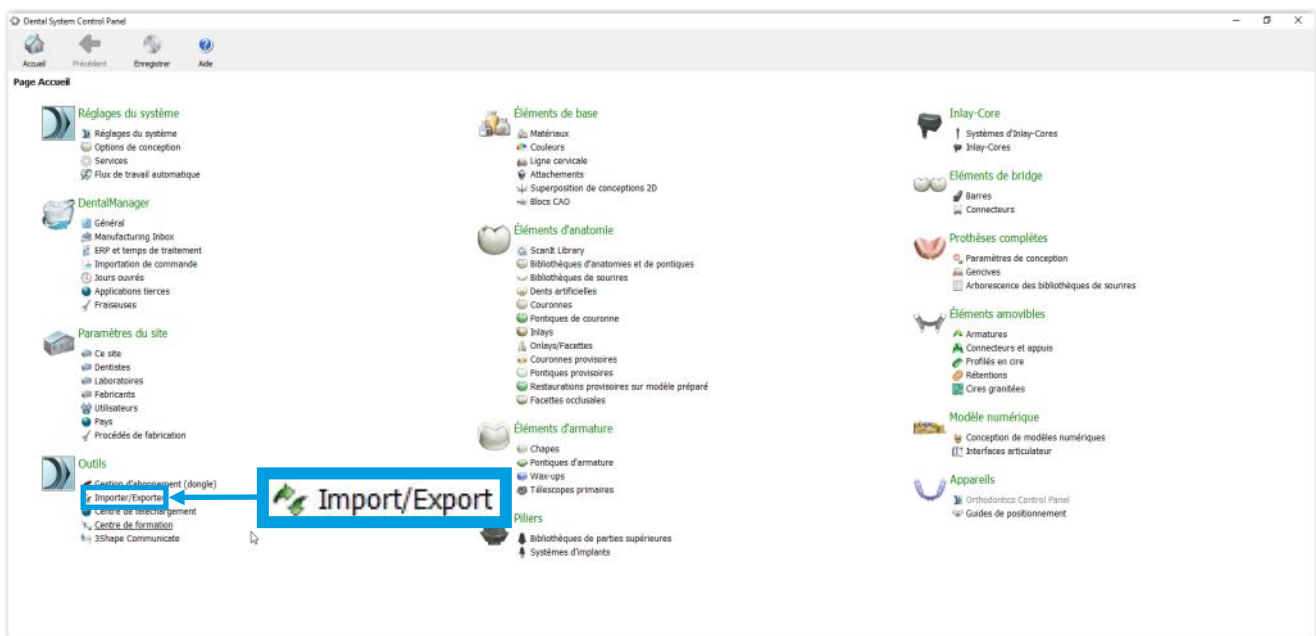

- Read the file: "Anthogyr-Labside-OPFLEX403\_503-SA-DT-PA-3Shape2018.dme" or  
"Anthogyr-Labside-OPFLEX403\_503-SA-DT-PA-3Shape2019.dme" or  
"Anthogyr-Labside-OPFLEX403\_503-SA-DT-PA-3Shape2020.dme".

- Then close the control panel.

- Restart the dental CAD system.

**Use:**

- Open an order form and select a tooth.
- Click on the abutment button on select :
  - the category : "Anthogyr-Labside-DT" or "Anthogyr-Labside-SA" (1)
  - the system : "Anthogyr Flexibase BL (DT)" or "Anthogyr Flexibase BL (SA)". (2)



- Then select the needed kit: after that start the scan (3):

