


Torq Control®

Universal Torque Wrench

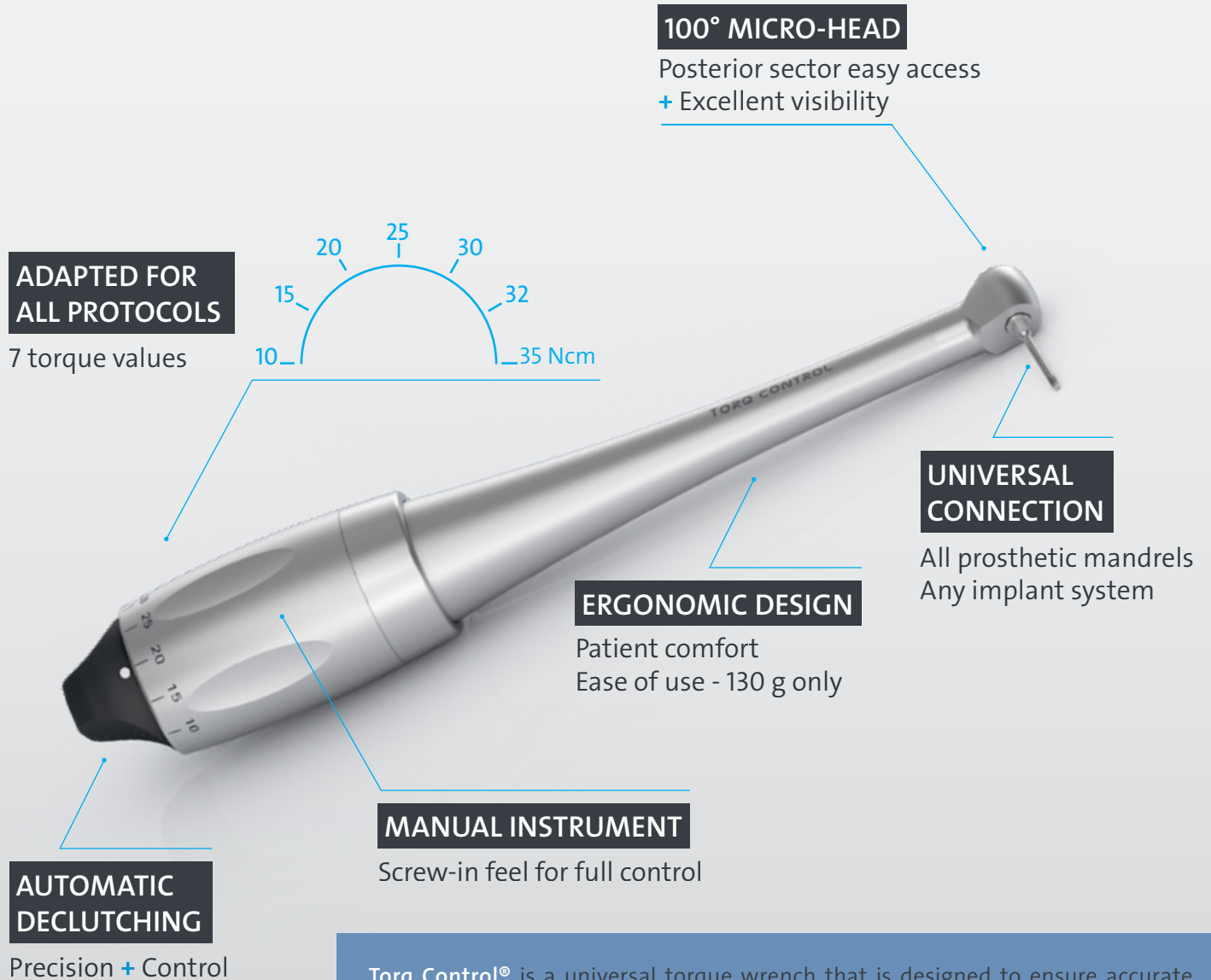

Torq Control® is a universal torque wrench that is designed to ensure accurate tightening regardless of the type of implant connection or accessibility challenges.

Tightening accuracy is crucial for successful implant treatment and Torq Control® has been specifically designed to meet these needs for all prosthetic procedures, with added control through automatic declutching.

Now available in North America! The original Torq Control is a must-have, especially for full-arch restorations.

Contact your local Anthogyr® representative for more information.

COMPATIBLE WITH ALL IMPLANT SYSTEMS!

- Universal connection for all dental implant systems
- 7 torque values, from 10 to 35 Ncm, adapted for all protocols

DESIGNED FOR EFFICIENCY

- Micro-head, 100° angulated for easy access and optimal visibility, even in posterior sectors

MAXIMUM COMFORT OF USE

- Manipulation 100% outside of mouth
- No lateral movement

SO EASY TO HANDLE

- Only 130 g for a better freedom of movement

CONTROL FEATURES


- High precision torque tightening
- Automatic declutching when the selected torque is reached

EASY MAINTENANCE

- One-piece design
- Sterilizable in autoclave, 135°

GUARANTEE OF QUALITY

- Made in the Mont-Blanc Valley, France, since 1997
- Top precision with more than 50 stages of high-quality manufacturing

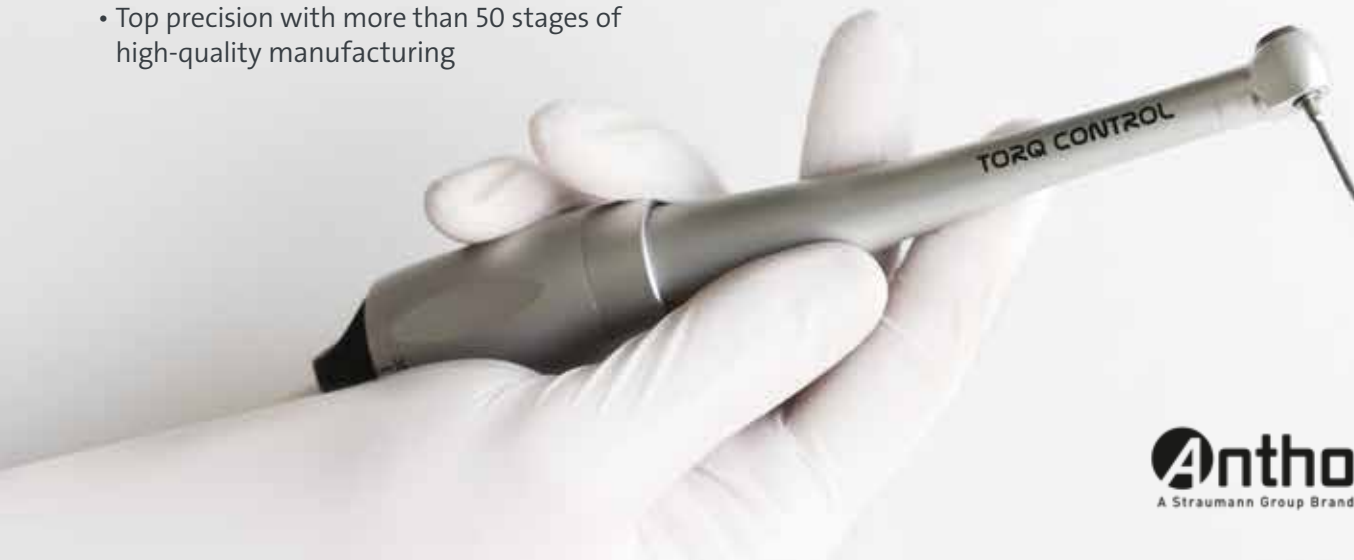

Anthogyr
A Straumann Group Brand

Straumann North American Headquarters

Straumann USA, LLC
60 Minuteman Road
Andover, MA 01810
Phone 800/448 8168 (US) • 800/363 4024 (CA)
Fax 978/747 2490
www.straumann.us • www.straumann.ca

ifu.anthogyr.com

© Straumann USA, LLC 2024. All rights reserved.

Straumann® and/or other trademarks and logos from Straumann® that are mentioned herein are the trademarks or registered trademarks of Straumann Holding AG and/or its affiliates. All rights reserved.