

INTEGRAL

guided surgery

SMILE IN A BOX® SERVICES

Focus on your passion - We take care of the rest



smileinabox.com

WALK ON A TRUSTED PATH

INTEGRAL

guided surgery



WALK ON A TRUSTED PATH

INTEGRAL is a fully guided surgery solution developed to address all your guided cases with the flexibility you need.

The **INTEGRAL** range gives you the freedom to match your clinical needs as closely as possible. Smart kits and instruments are designed to enhance the accuracy and ergonomics of your surgical

procedures. **INTEGRAL** is an open solution that allows you to integrate major workflows.

Enjoy the predictability of guided surgery. Choose **INTEGRAL** and its distinctive design developed to increase your surgical efficiency.

USLIT.3010 V1 PMR 2/24 ifu.anthogyr.com

Photos credits: Anthogyr, Unsplash. All rights reserved, actual products may vary.


Anthogyr
A Straumann Group Brand

www.anthogyr.us

Anthogyr
A Straumann Group Brand

Anthogyr
A Straumann Group Brand

Anthogyr | INTEGRAL guided surgery - Walk on a trusted path

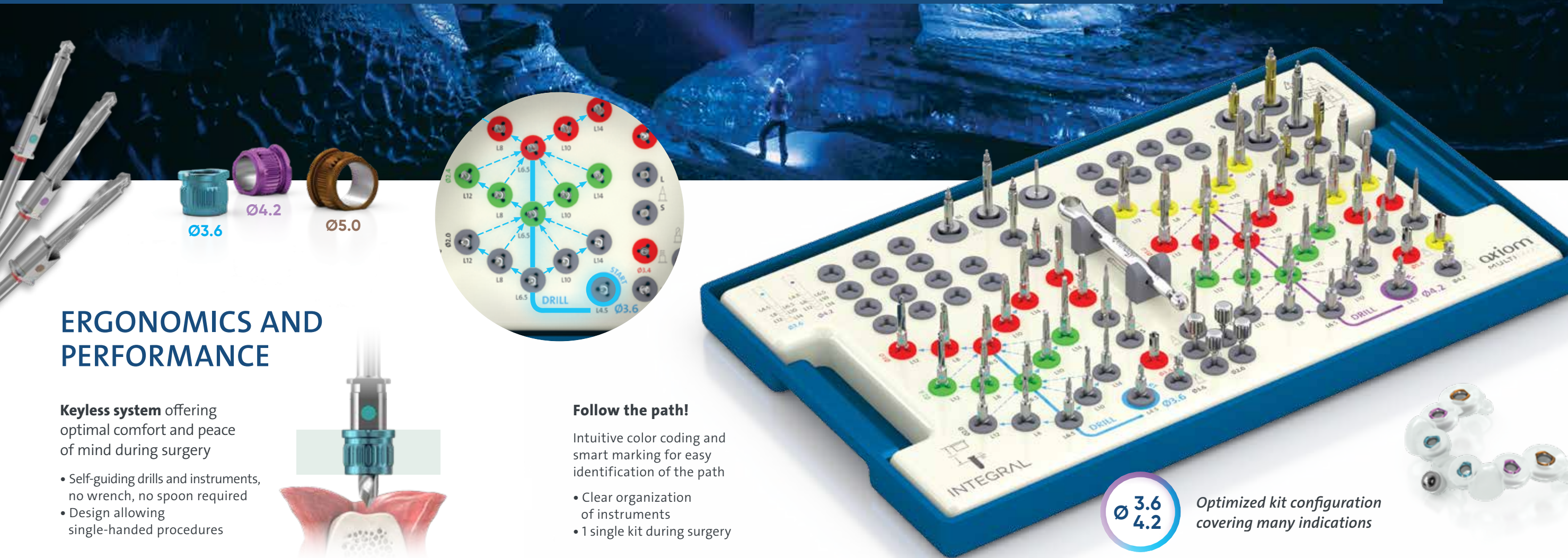

ERGONOMICS AND PERFORMANCE

Keyless system offering optimal comfort and peace of mind during surgery



- Self-guiding drills and instruments, no wrench, no spoon required
- Design allowing single-handed procedures



Follow the path!


Intuitive color coding and smart marking for easy identification of the path

- Clear organization of instruments
- 1 single kit during surgery

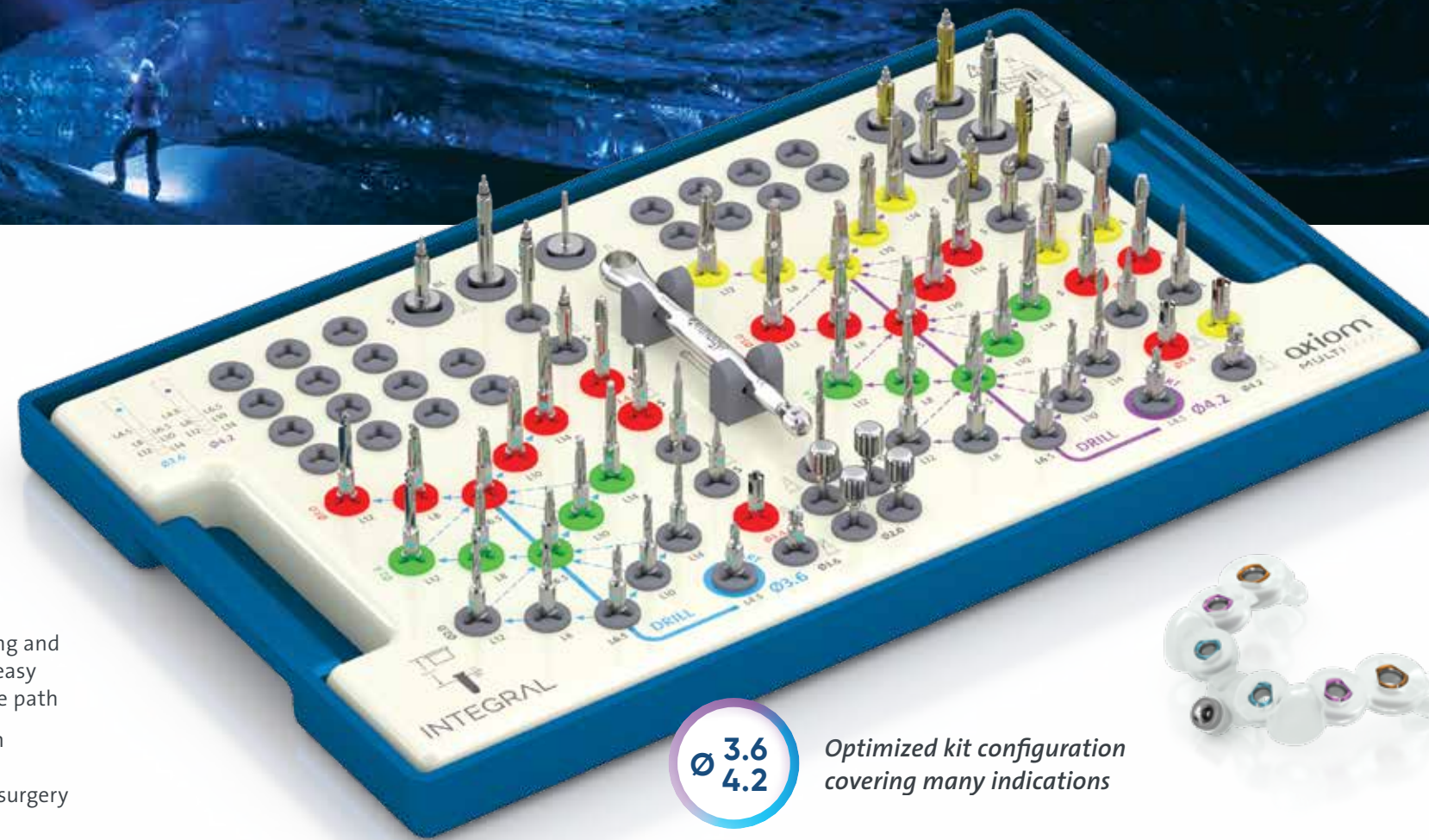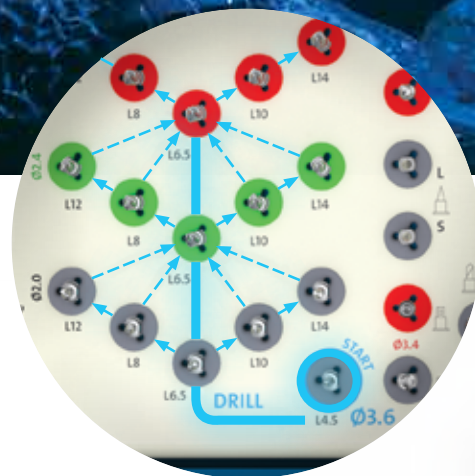

Ø 3.6
Ø 4.2

Optimized kit configuration covering many indications

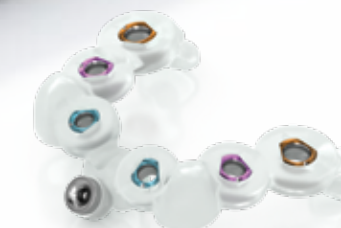

FREEDOM

Open compatibility with major planning software


- Seamless integration into dental offices and optimized communication with dental labs



ACCURACY

V-shape sleeves with implant indexing

- Optimized design for interdental spaces
- Unique placement of the sleeves in the surgical guide
- Markers for appropriate final implant placement in depth and axis
- Designed to provide comfortable handling

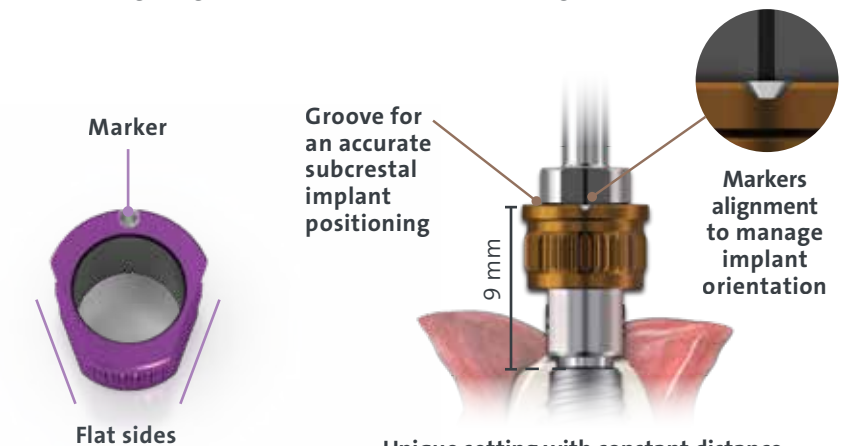


FLEXIBILITY

Enjoy the **flexibility of an exhaustive range of 4 kits**, to cover main clinical indications with Axiom® implants.

axiom®

- Axiom X3® and Axiom® PX
- Ø3.4 mm to Ø4.6 mm
- 6.5 mm to 14 mm lengths

Ø 3.6
Ø 4.2

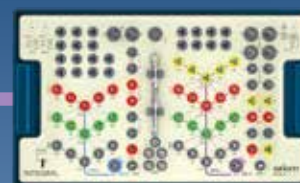


- **Main clinical situations** ♥
- Axiom® Ø3.4 / 4.0 mm

Ø3.6




- **Narrow interdental spaces**
- Axiom® Ø3.4 mm

Ø4.2

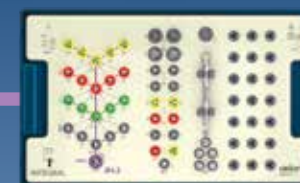

- **Intermediate interdental spaces**
- Axiom® Ø3.4 / 4.0 mm

Ø5.0

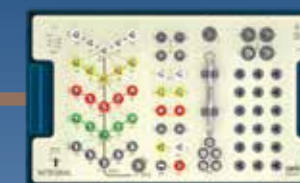

- **Large interdental spaces**
- Axiom® Ø3.4 / 4.0 / 4.6 mm