


The premier aligner

Engineered for precision.

A unique combination engineered to move teeth more predictably.

What does it take to move teeth? It's simple, really: consistent, continuous force. ClearCorrect™ aligners are uniquely designed to deliver exactly that, because only ClearCorrect aligners feature a flat trimline above the gingival margin, combined with our proprietary ClearQuartz™ tri-layer aligner material and new customizable features. The resulting aligner is remarkably retentive and effective, and doesn't sacrifice comfort, clarity, or resilience as it helps you achieve optimal results in a wide variety of cases.

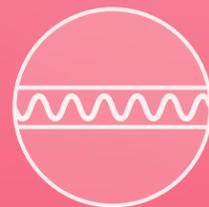

ClearQuartz
Tri-layer material




High, flat trimline



Customizable features

ClearQuartz™ tri-layer material

ClearQuartz features an elastomeric layer between two resilient, low-porosity materials. The outer layers allow the aligner to be tough and stain-resistant, while the inner layer allows the aligner to apply gentle, consistent, continuous force.

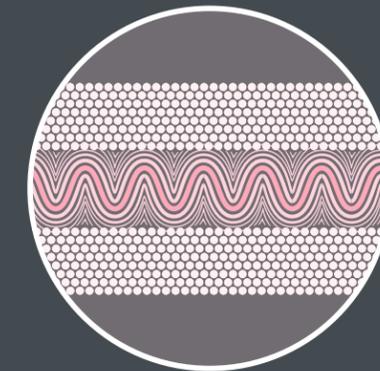


Low-porosity material
For toughness and stain-resistance

Elastomeric layer
Applies gentle, consistent, and continuous force

ClearQuartz retains 10x more of its initial force than competitors¹

and retains its shape throughout the wear time of the aligner².

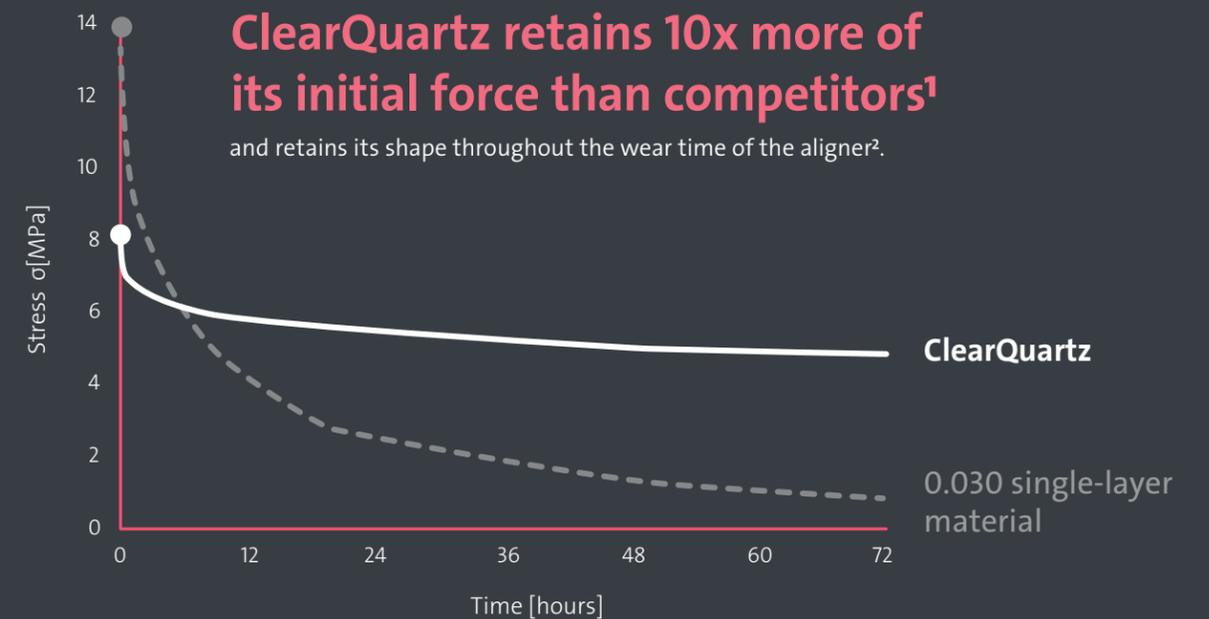


¹ Data on file. Over 72 hours.

² Data on file vs 0.30 single layer material.

The ClearCorrect trimline: Proven to optimize fit and retention.


With a smooth, unscalloped trimline that extends past the gingival margin, ClearCorrect aligners are more than 2.5x more retentive than scalloped aligners, and are engineered to achieve more complex tooth movements, even while reducing the need for engagers.

ClearCorrect aligners are over 2.5x more retentive vs. scalloped trimline³

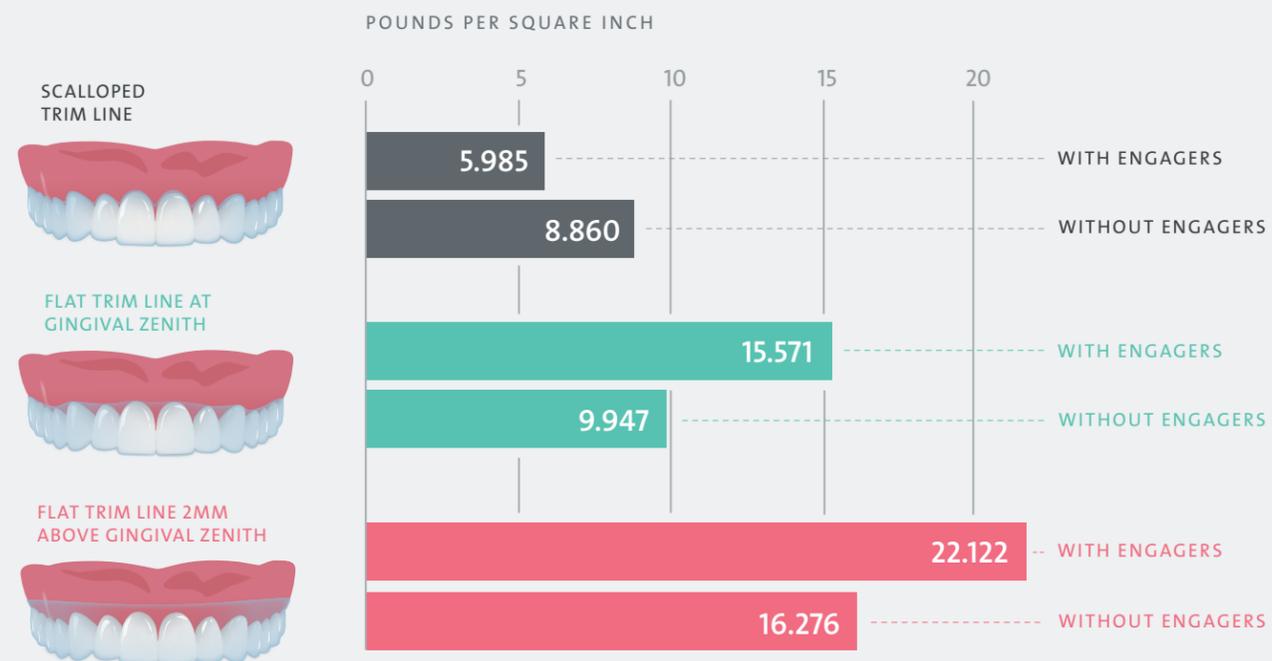

Customizable to suit your patient's needs.

All patients are different, and each aligner is custom-made to fit, but with specialized options for engagers and cutouts, ClearCorrect is giving you even more ways to tailor your patient's treatment precisely to their needs.

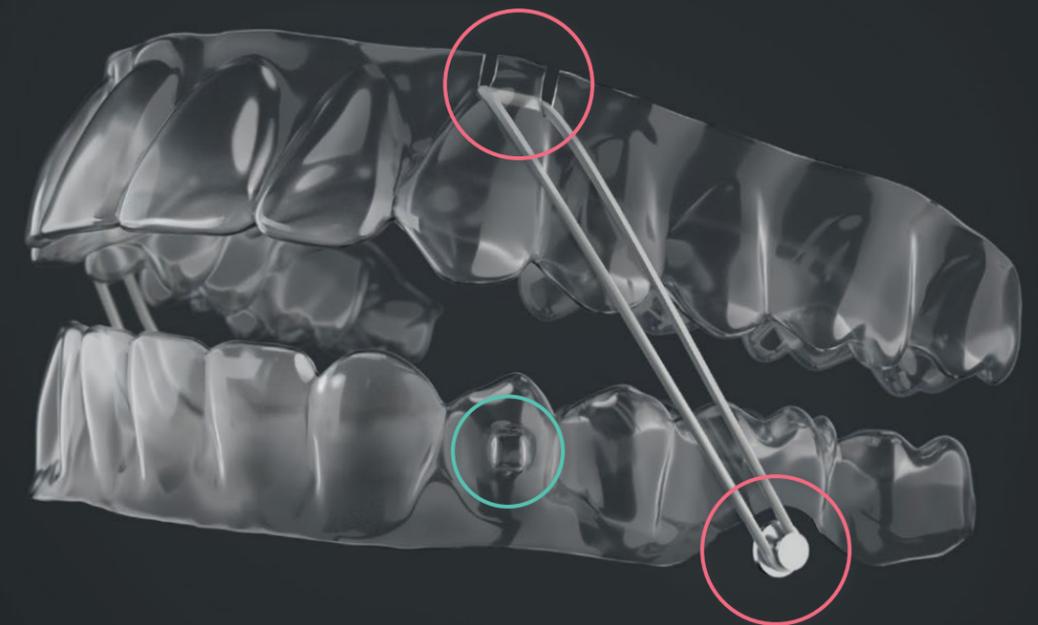

ENGAGERS

Engager sizes are now available in 2mm, 3mm, and 4mm. Decide what size and orientation are best for your case.


CUTOUTS AND ELASTICS

With cutouts for elastics ClearCorrect offers you more control to perform more complex movements. Confidently treat teens, Class II, and Class III malocclusions.

- Select a double vertical slit or half-moon button cutout
- Prescribe them in the Doctor Portal and approve them in ClearPilot™
- Receive your aligners with pre-made cutouts for your patient to start wearing according to your treatment plan.

For the ultimate patient experience.

A more flexible aligner that adapts to teeth better, with **33% less initial force**⁴ and remarkable clarity and stain resistance? No wonder patients like our aligners so much.

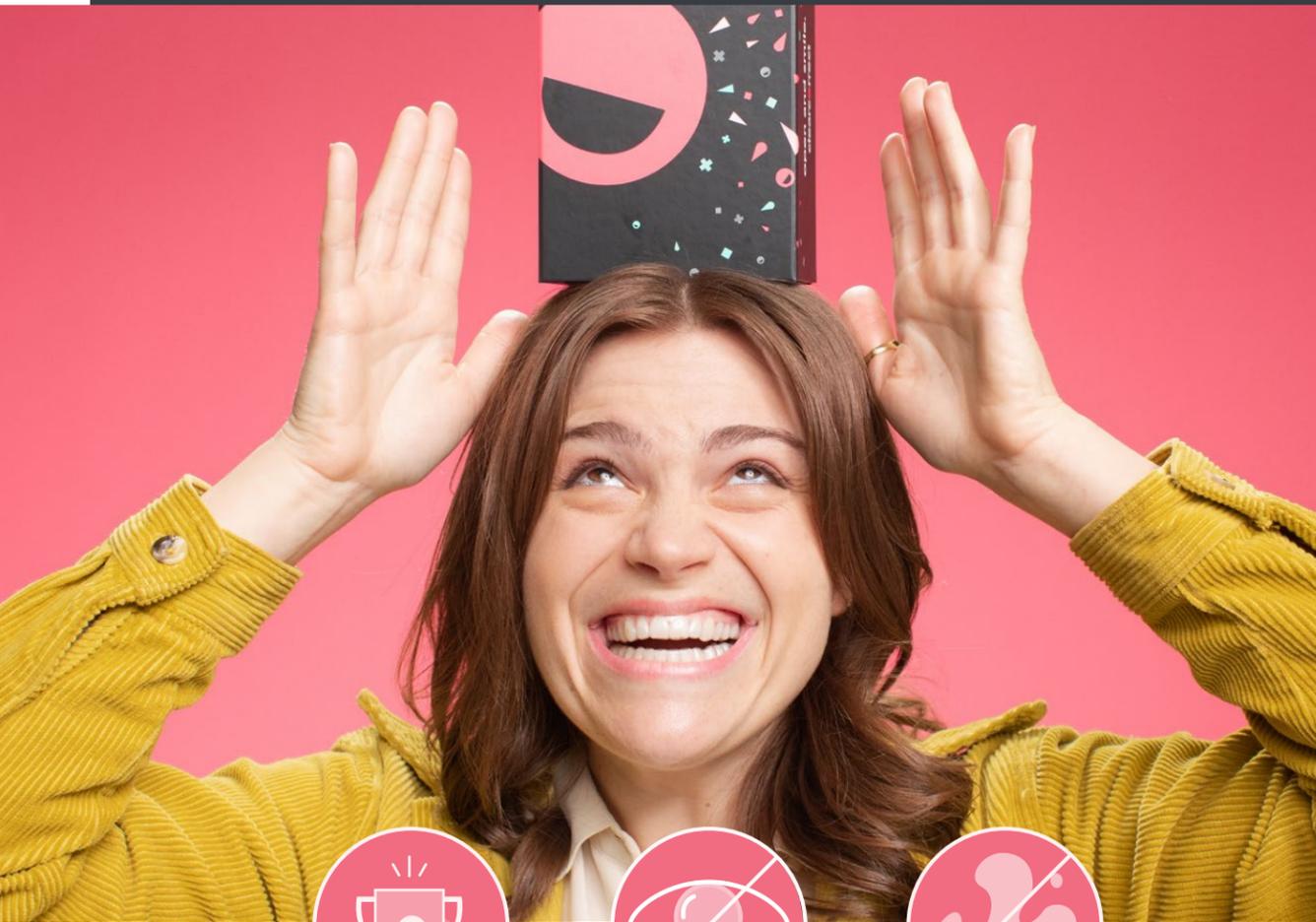

Improved comfort
1/3 less initial force



Nearly invisible



Stain resistant

Don't take it from us. Take it from them.



Dr. Rebecca Komischke

ClearQuartz has made aligner treatment so much more comfortable for my patients. Some were in mid-treatment when the new material was launched. They were amazed by the new feel and fit.

Dr. Ly T. Mar

It's been beautiful...patients are happy, they love it. A lot less emergencies, a lot more efficient product: better grip on the teeth, better movement, more predictable movement...we are doing a lot of ClearCorrect products right now.



The aligners are only part of the story.

As impressive as they are, our aligners are only part of what makes ClearCorrect your partner in orthodontics. The ClearCorrect system combines our premier aligner with an empowering digital workflow and comprehensive educational tools to deliver the support you need, with the freedom you want and the products that you and your patients deserve.





Register & Submit.

dr.clearcorrect.com

Learn more.

learn.clearcorrect.com

Get support.

support@clearcorrect.com

See real cases.

results.clearcorrect.com

clearcorrect.com



clearcorrect
A Straumann Group Brand