



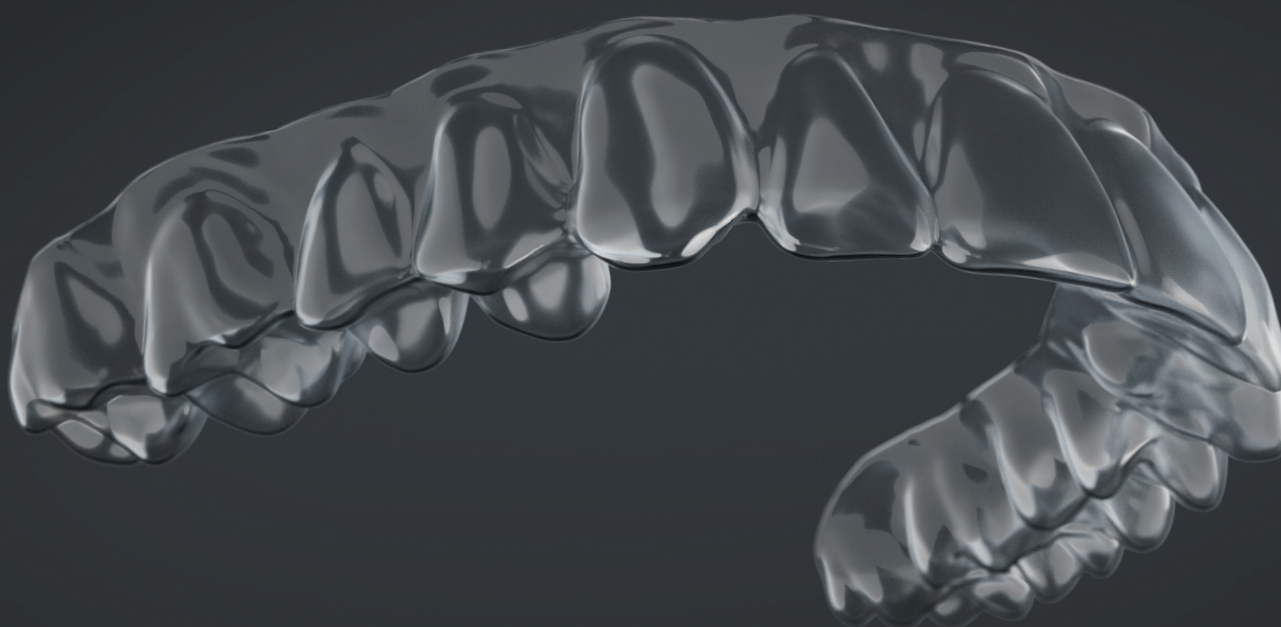


The premier aligner


Engineered for precision.

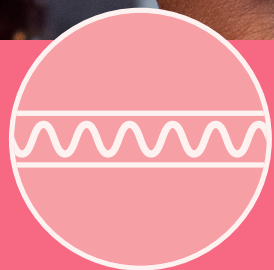

ClearQuartz
Tri-layer material

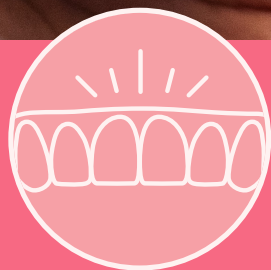

High, flat trimline

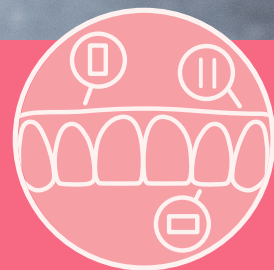

Customizable features


A unique combination engineered to move teeth more predictably.

What does it take to move teeth? It's simple, really: consistent, continuous force. ClearCorrect® aligners are uniquely designed to deliver exactly that, because only ClearCorrect aligners feature a flat trimline above the gingival margin, combined with our proprietary ClearQuartz™ tri-layer aligner material and new customizable features. The resulting aligner is remarkably retentive and effective, and doesn't sacrifice comfort, clarity, or resilience as it helps you achieve optimal results in a wide variety of cases.




ClearQuartz™ tri-layer material

ClearQuartz features an elastomeric layer between two resilient, low-porosity materials. The outer layers allow the aligner to be tough and stain-resistant, while the inner layer allows the aligner to apply gentle, consistent, continuous force.

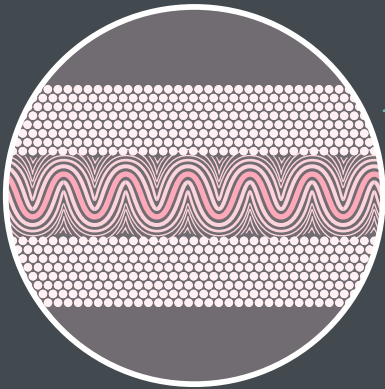

Low-porosity material
For toughness and stain-resistance

Elastomeric layer
Applies gentle, consistent, and continuous force

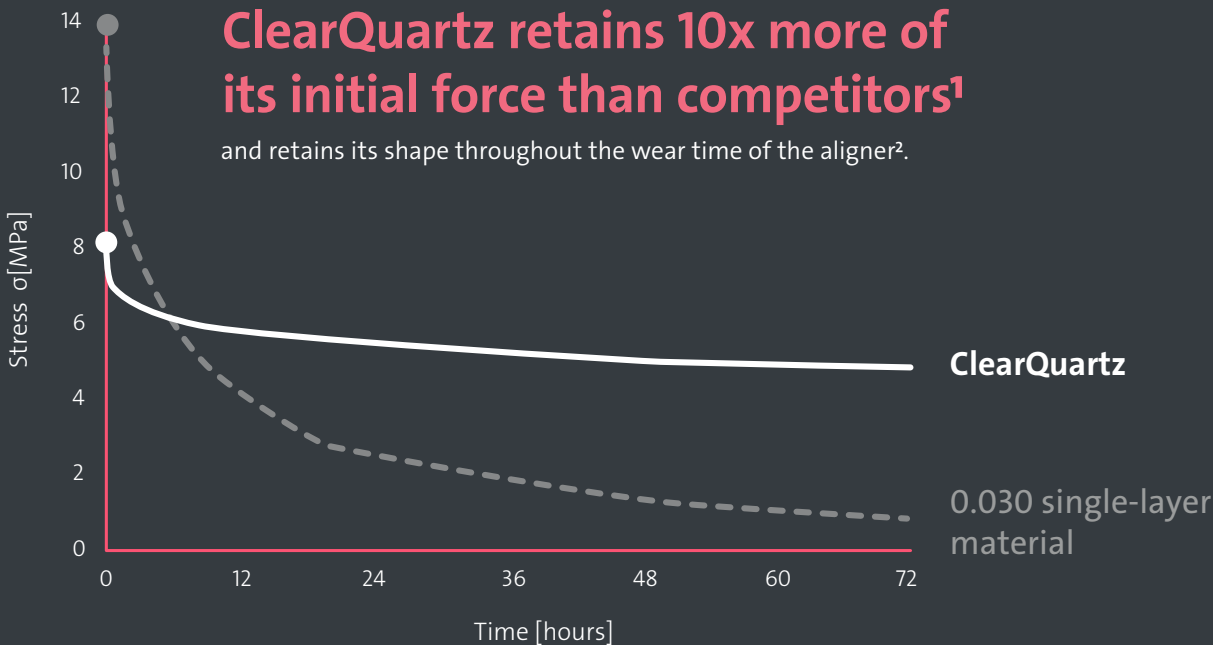

ClearQuartz retains 10x more of its initial force than competitors¹
and retains its shape throughout the wear time of the aligner².

¹ Data on file. Over 72 hours.
² Data on file.

Stain-resistant material

ClearCorrect has been demonstrated to be stain-resistant versus the leading aligner brand. In a simulated test, ClearCorrect remains clearer with 4 different staining agents: Cola, Red Wine, Coffee, and Mustard³. This result may support the hypotheses that the elastomeric layer (middle layer in ClearCorrect's aligner) is responsible for the absorption of staining compounds, as suggested by previous studies.^{4,5}



COLA

WINE

ClearCorrect

ClearCorrect

Invisalign®

Invisalign®

MUSTARD

COFFEE

ClearCorrect

ClearCorrect

Invisalign®


Invisalign®

³ Data on file.
⁴ Liu, C.-L. et al. Colour stabilities of three types of orthodontic clear aligners exposed to staining agents. Int. J. Oral Sci. 8, 246–253 (2016).
⁵ Bernard G, Rompré P, Tavares JR, Montpetit A. Colorimetric and spectrophotometric measurements of orthodontic thermoplastic aligners exposed to various staining sources and cleaning methods. Head Face Med. 16, 1–11 (2020).
Invisalign is a registered trademark of Align Technology.

The ClearCorrect trimline: Proven to optimize fit and retention.

With a smooth, unscalloped trimline that extends past the gingival margin, ClearCorrect aligners are more than 2.5x more retentive than scalloped aligners, and are engineered to achieve more complex tooth movements, even while reducing the need for engagers.

ClearCorrect aligners are over 2.5x more retentive vs. scalloped trimline⁶

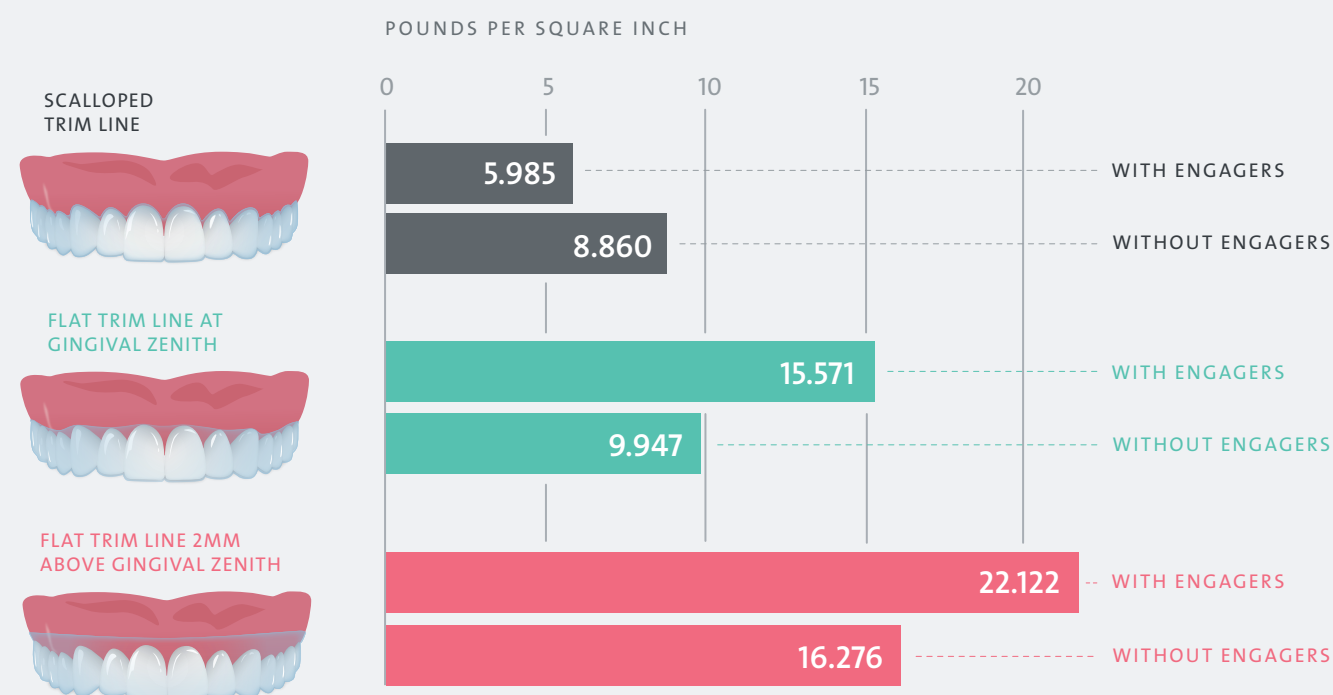

For the ultimate patient experience.

A more flexible aligner that adapts to teeth better, with **33% less initial force⁷** and remarkable clarity and stain resistance? No wonder patients like our aligners so much.



Improved comfort
1/3 less initial force




Nearly invisible



Stain resistant

Customizable to suit your patient’s needs.

All patients are different, and each aligner is custom-made to fit, but with specialized options including engagers, cutouts for elastics, and bite ramps, ClearCorrect is giving you even more ways to tailor your patient’s treatment precisely to their needs.

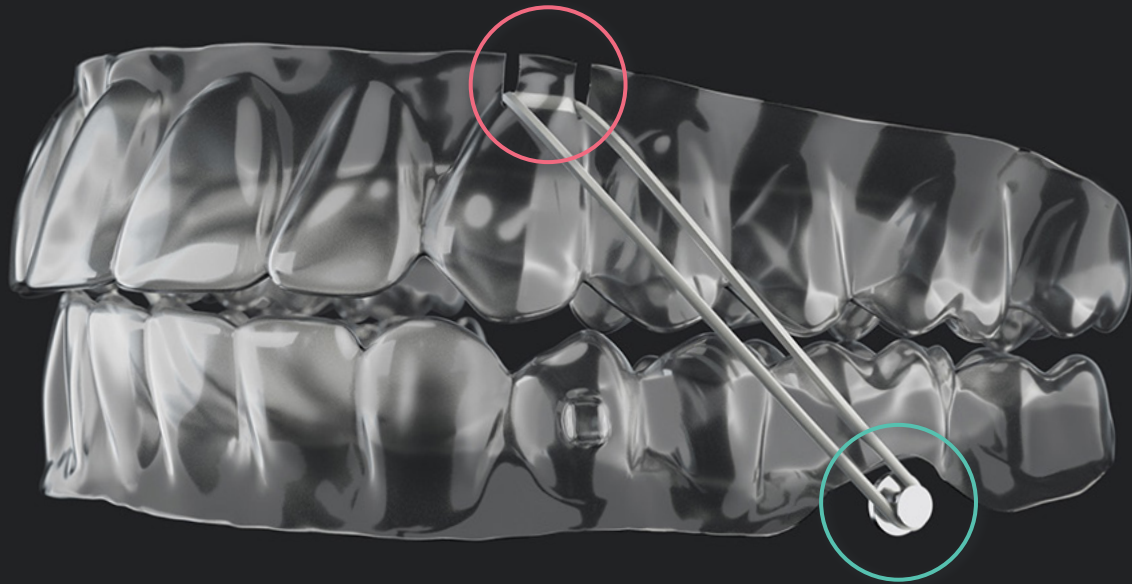

CUTOUTS FOR ELASTICS

With cutouts for elastics ClearCorrect offers you more control to perform more complex movements. Confidently treat teens, Class II, and Class III malocclusions.

BUTTON CUTOUTS

- Half-moon design
- Fits buttons up to 6mm in diameter
- Elastic anchored directly to the button bonded on the tooth
- Available for buccal placement on any teeth, or lingual placement on molars or premolars
- Staged from first or any preferred odd step

SLIT CUTOUTS

Double vertical slit cutout design engineered to optimize treatment with elastics in combination with ClearCorrect’s tri-layer material and high, flat trimline.

- Designed for elastic to attach directly to the aligner
- Available for buccal or lingual placement on any teeth
- Staged from first or any preferred odd step

ClearCorrect offers customized engagers, cutouts, and bite ramps—manufactured ready for wear according to your patient’s unique treatment plan. You can request these features with any treatment option, independent of—and in tandem with—the treatment preferences you set in the Doctor Portal.



1. Add your preferences




2. Submit your prescription



3. Review, request edits, or approve your setup



4. Receive the aligners ready to wear

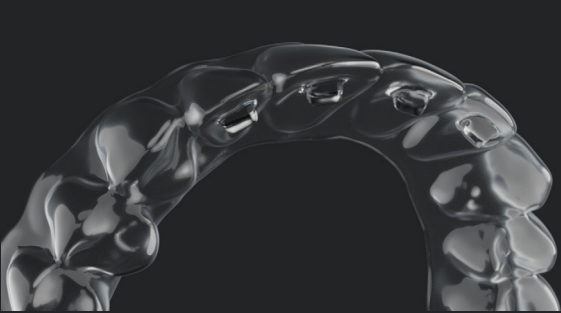

BITE RAMPS

Expand your treatable cases to include Deep Bites and Crossbites with bite ramps from ClearCorrect.

- Facilitates intrusion of lower anteriors
- Disocclude the posterior teeth to assist bite correction
- Customizable to your patient needs in 2 different sizes: 3mm and 5mm depth
- Available in the configuration of your preference: 2x2, 3x3, or only 3s
- Bite ramps planned to keep contact with lower anteriors



ENGAGERS

With more choices to select from, customize your patient’s treatment with an extended selection of engagers. Decide what size and orientation are best for your case.

- All engagers are rectangular in shape
- Available in horizontal and vertical orientation
- Available in 2, 3 and 4mm sizes
- Available in 0.75*, 1.00, and 1.25** mm depths



The aligners are only part of the story.

As impressive as they are, our aligners are only part of what makes ClearCorrect your partner in orthodontics. The ClearCorrect system combines our premier aligner with an empowering digital workflow and comprehensive educational tools to deliver the support you need, with the freedom you want and the products that you and your patients deserve.

Don't take it from us. Take it from them.



Dr. Rebecca Komischke

ClearQuartz has made aligner treatment so much more comfortable for my patients. Some were in mid-treatment when the new material was launched. They were amazed by the new feel and fit.




Dr. Ly T. Mar

It's been beautiful...patients are happy, they love it. A lot less emergencies, a lot more efficient product: better grip on the teeth, better movement, more predictable movement...we are doing a lot of ClearCorrect products right now.



Dr. Daniel Neves

The release of ClearQuartz reaffirms the fact that the Straumann Group is committed in bringing the maximum of innovation and technology to the orthodontic field.

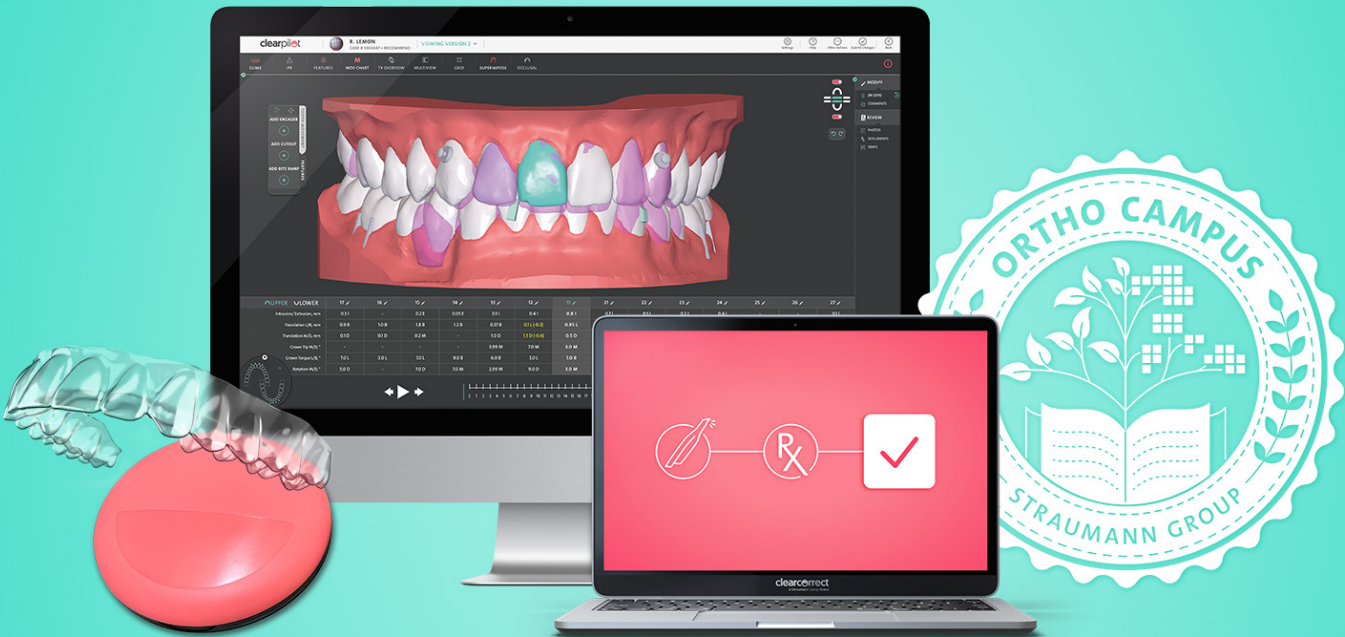



Register & Submit.

dr.clearcorrect.com

Learn more.

learn.clearcorrect.com

Get support.

support@clearcorrect.com

See real cases.

results.clearcorrect.com

clearcorrect.com



© ClearCorrect 2022 | 480.530_en_06 07/05 | Approved for global use.

Not all features are available in all markets.

clearcorrect
A Straumann Group Brand