


Customizable to suit your patient's needs.

All patients are different, and each aligner is custom-made to fit, but with specialized options including engagers, cutouts for elastics, and bite ramps, ClearCorrect® is giving you even more ways to tailor your patient's treatment precisely to their needs.

CUTOUTS FOR ELASTICS

With cutouts for elastics ClearCorrect offers you more control to perform more complex movements. **Confidently treat teens, Class II, and Class III malocclusions.**

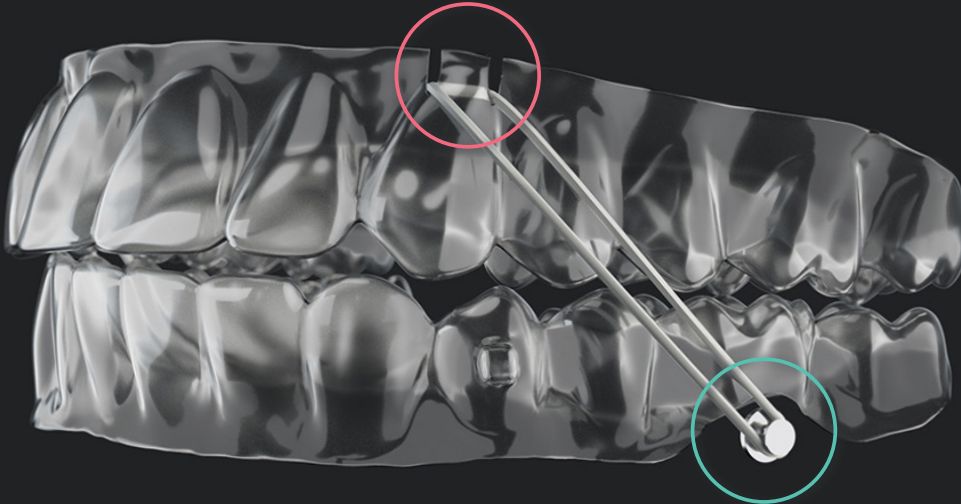


BUTTON CUTOUTS

- Half-moon design
- Fits buttons up to 6mm in diameter
- Elastic anchored directly to the button bonded on the tooth
- Available for buccal placement on any teeth, or lingual placement on molars or premolars
- Staged from first or any preferred odd step

SLIT CUTOUTS

Double vertical slit cutout design engineered to optimize treatment with elastics in combination with ClearCorrect's tri-layer material and high, flat trimline.


- Designed for elastic to attach directly to the aligner
- Available for buccal or lingual placement on any teeth
- Staged from first or any preferred odd step

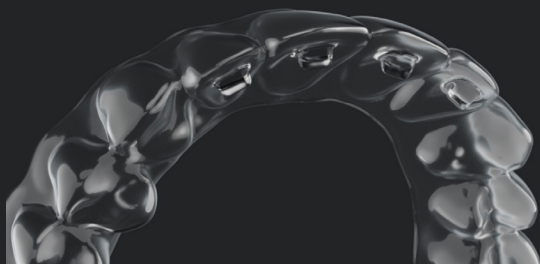

BITE RAMPS

Expand your treatable cases to include Deep Bites and Crossbites with bite ramps from ClearCorrect.

- Facilitates intrusion of lower anteriors
- Disocclude the posterior teeth to assist bite correction
- Customizable to your patient needs in 2 different sizes: 3mm and 5mm depth
- Available in the configuration of your preference: 2x2, 3x3, or only 3s
- Bite ramps planned to keep contact with lower anteriors

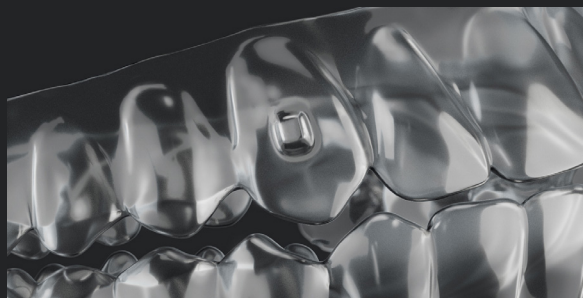

ENGAGERS

With more choices to select from, customize your patient's treatment with an extended selection of engagers. Decide what size and orientation are best for your case.

- All engagers are rectangular in shape
- Available in horizontal and vertical orientation
- Available in 2, 3 and 4mm sizes
- Available in 0.75*, 1.00, and 1.25** mm depths

*Only on horizontal engagers

**Only on vertical engagers

ClearCorrect offers customized engagers, cutouts, and bite ramps—manufactured ready for wear according to your patient's unique treatment plan. You can request these features with any treatment option, independent of—and in tandem with—the treatment preferences you set in the Doctor Portal.



1. Add your preferences



2. Submit your prescription



3. Review, request edits, or approve your setup



4. Receive the aligners ready to wear



Scan to learn
more today!

480.990_en_02 | Approved for global use.

Not all features are available in all markets.

clearcorrect
A Straumann Group Brand