

Introducing ICam

True Photogrammetry



Revolutionizing Implantology through Photogrammetry

The ICam is a handheld “camera unit” consisting of four cameras and one projector. ICam combines photogrammetric and structured light scanning techniques to capture 3D data. Using the ICamBodies, high precision mechanical parts with a proprietary target arrangement, helps determine the position and orientation of implants.



INNOVATIVE TECHNOLOGY

ICam is equipped with 4 moving cameras, compared to other stereo-photogrammetry or 2 stationary cameras available on the market



STREAMLINED WORKFLOW

In combination with an IOS like Straumann SIRIOS™, ICam can reduce the need for physical dental models, making your process more efficient



PEACE OF MIND

Easy to learn, predictable. Robust limited manufacturer warranty and support options keep you in production

CREATE YOUR SET OF ICAMBODIES


Set of 8 or 16 available.

BUILT TO LAST

Made from titanium to withstand torquing and autoclaving.

CAPTURE DIGITALLY

In lieu of a traditional implant impression, insertion of ICamBodies and the capture of gingiva are rapid and support patient comfort.

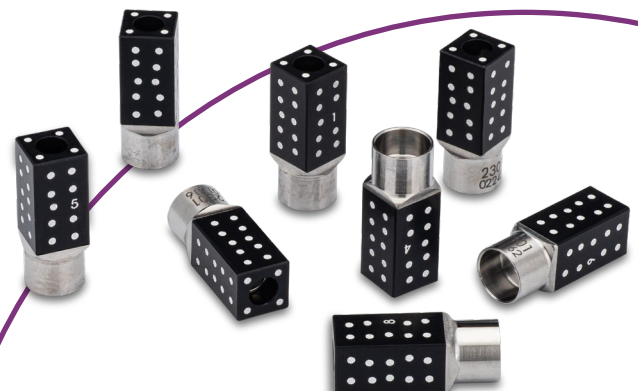


Enhance your Digital Journey

Discover the power of the Straumann digital solution, a comprehensive portfolio designed to grow your practice with full arch treatment. This all-encompassing suite includes advanced IOS technology, cutting-edge photogrammetry, 3D printers, milling machines, and a wide range of discs and resins. The Straumann digital solution provides everything you need to seamlessly integrate and scale full arch procedures, supporting both clinical efficiency and patient outcomes.

Digital Workflow* End-to-End

See how a streamlined workflow can help optimize treatment processes for fixed full arch. Plus, an In lab validated custom abutment workflow, as well as comprehensive Straumann Digital Solutions including planning, design, milling service with UN!Q, supporting precision and reliability from start to finish.

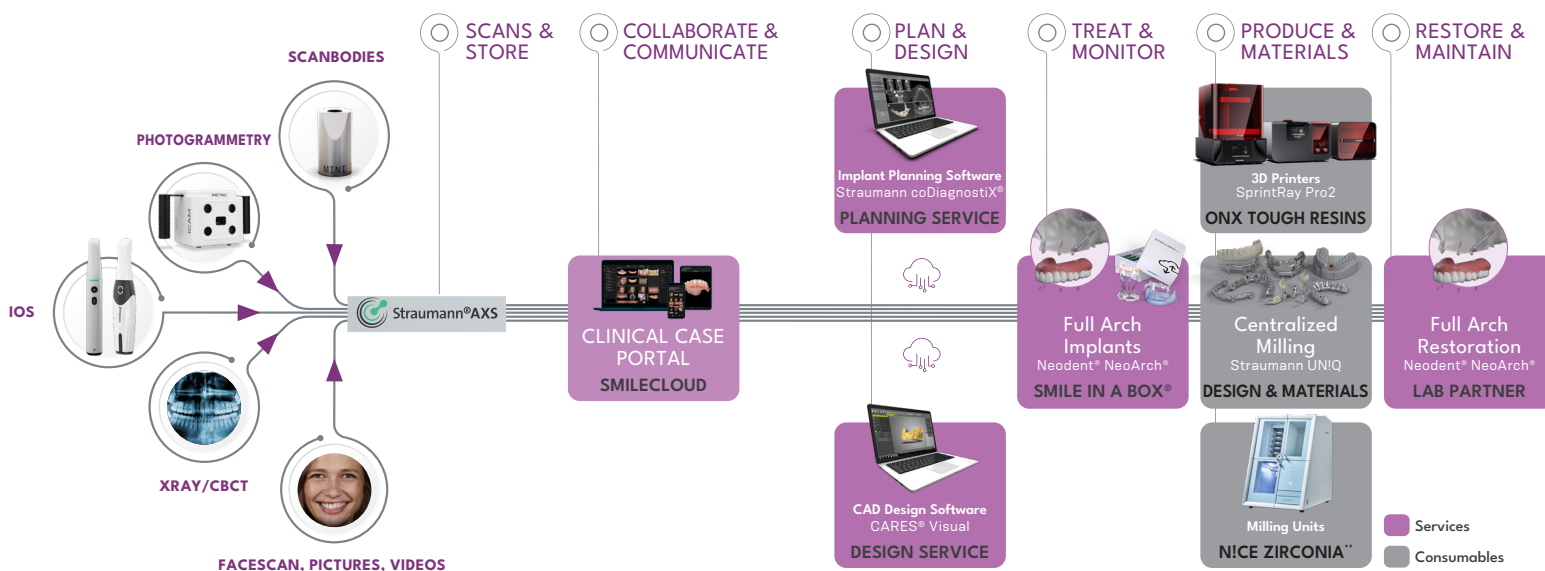

*All digitally designed copings and/or crowns to be used with the Neodent Titanium Base Abutment System are intended to be sent to Straumann for manufacture at a validated milling center.

**Not all products depicted are available in all markets.