

The Straumann® Bone Level Tapered Implant Ø2.9 mm



One system,
one final abutment torque

The Bone Level Tapered Ø2.9 mm
is indicated for the following:

- maxillary lateral incisors
- mandibular incisors

Roxolid® is a material with
excellent mechanical properties

Apically tapered implant body
is designed to allow
underpreparation to support
primary stability in soft bone

Variobase® for crown
for single-unit restorations
of Straumann SC Implants

Small CrossFit® connection
provides simplified accuracy
and provides confidence for
component positioning

SLActive® surface allows
faster osseointegration* and
increased predictability




*Compared to SLA

More than a small diameter.

Slim yet strong.

Over **6 million people**, approximately 2% of the population, in the United States have one or both of the maxillary lateral incisors missing.¹ The reduced interdental space of the maxillary and mandibular incisors can be challenging to treat. The Straumann® BLT Ø2.9 mm implant is designed to treat these patients while providing peace of mind.

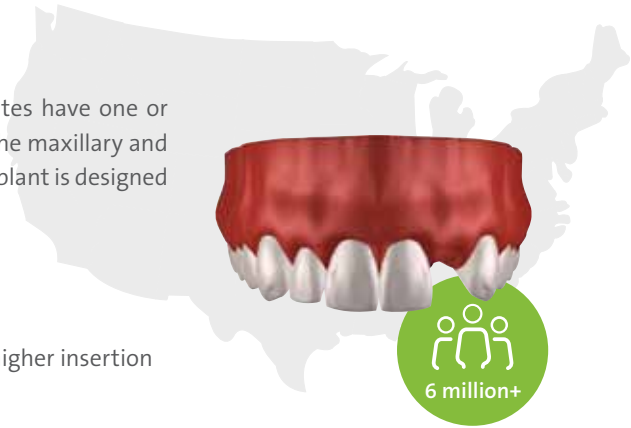


HIGH IMPLANT-ABUTMENT STRENGTH

The Bone Level Tapered Ø2.9 mm implant shows high mechanical strength and higher insertion torque for final abutments.

| | STMN BLT Ø2.9 mm | Nobel Active Ø3.0 mm Nobel Biocare | Osseospeed Ø3.0 mm Astra Tech |
|--------------|------------------|---------------------------------------|----------------------------------|
| Torque [Ncm] | 35 | 15* | 15** |

STRAUMANN® BONE LEVEL TAPERED Ø2.9 MM IMPLANT CLINICAL CASE

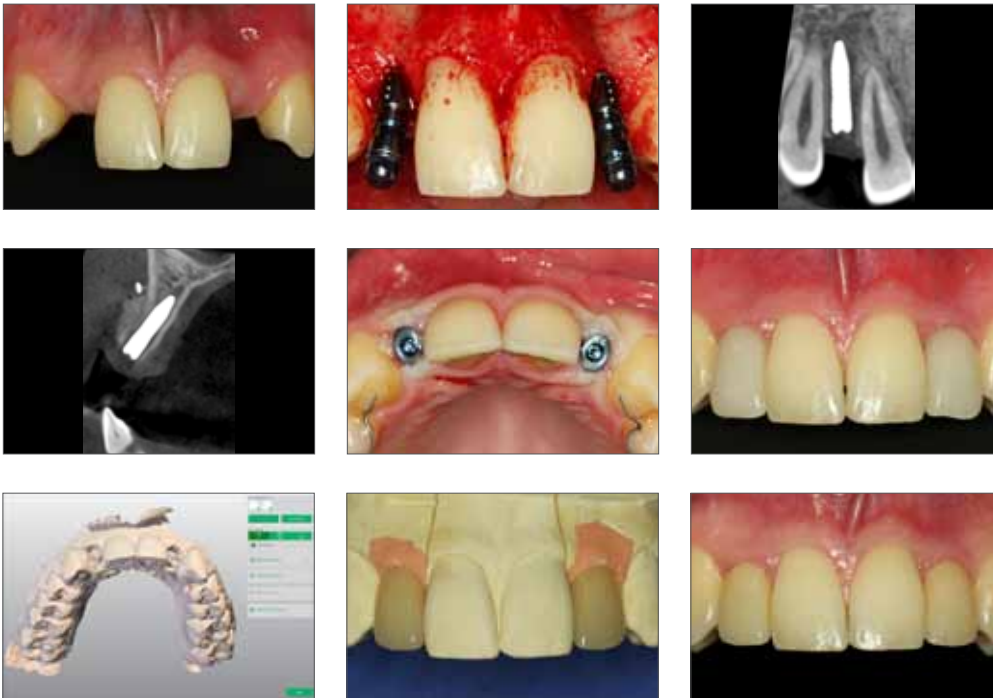

Case courtesy: Dr. med. dent. Kay Vietor

Data on file @Straumann
Test set up according to Norm ISO 14 801

¹Symons AL, Stritzel F, Stauation J. Anomalties associated with hypodontia of the permanent lateral incisor ans second premolar. J Clin Pediatr Dent. 1993 Winter; 17(2): 109-11.

*Prosthetic Torque Guide. <https://www.nobelbiocare.com/content/dam/Migration%20Assets/Documents/Discover/customer-service/Prosthetic%20Torque%20Guide.pdf>. Accessed December 20, 2016.

**Astra Tech Implant System Product Catalog Astra Tech Implant System for OsseoSpeed TX. http://www.dentsplyimplants.com/~/_/media/M3%20Media/DENTSPLY%20IMPLANTS/Product/1207012%20Product%20catalog.ashx?filetype=.pdf. Accessed December 20, 2016.