

STRAUMANN® SURGICAL MOTOR **PRO 2**

Designed to improve operability,
safety, wireless extensibility.*



* vs Straumann® Surgical Motor Pro

Product overview

STRAUMANN® SURGICAL MOTOR PRO 2

The Straumann® Surgical Motor Pro has evolved into Straumann® Surgical Motor Pro 2. The main focus was on further improving operability, safety and wireless extensibility.



Strong

- A robust motor with maximum torque value of 80 Ncm, suitable for all Straumann® implant types
- 2 year warranty on control unit and micromotor



Simple

- Large display with all information on one screen
- Natural daylight LED illumination
- Memory for 8 implant protocols
- User-friendly foot control allowing operation of all functions within the preset parameters without touching the control panel




SAFE

- “Torque limiter” stops the motor when the working torque reaches the preset torque value
- “Torque calibration” ensures accurate torque output for any contra-angle handpiece

OPERABILITY

There are a range of upgrades to ensure great comfort during implant treatment. A compact, shortened motor, (reduced by more than 2 cm) to improve ergonomics.

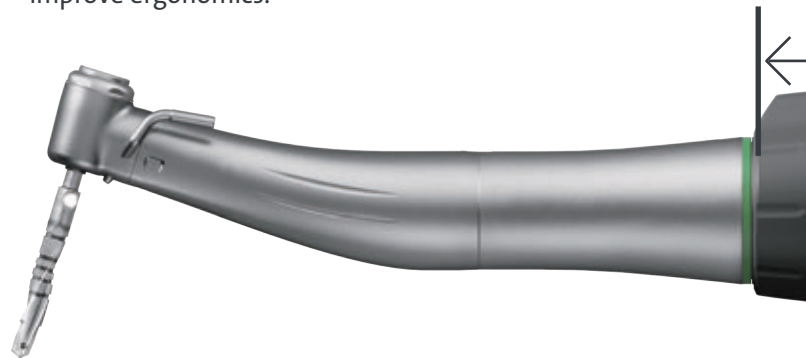

- 1 Cooling solution flow volume button
- 2 Program selection button
- 3 Forward/reverse button
- 4 Footpedal button





SAFETY



The combined usage of Straumann® Surgical Motor Pro 2 and the osseointegration monitoring device Straumann® Osseo 100+*, allows to measure contactless Implant Stability Quotient (ISQ). The 'Advanced Handpiece Calibration' (AHC) function is designed to provide accurate torque values.



EXTENSIVITY





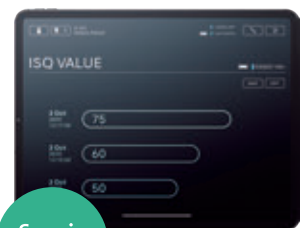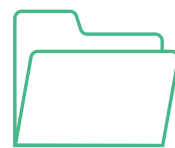


Enables the communication between different devices by using the wireless connection.

The osseointegration monitoring device, Straumann® Osseo 100+, can be connected to the Straumann Surgical Motor Pro 2 using Bluetooth®. This allows to transfer automatically the ISQ value.



* Straumann® Osseo 100+ might not be available in all the countries

STRAUMANN® SURGICAL MOTOR PRO 2

	Article Name	Article Number
	Straumann® Surgical Motor Pro 2 Contains: <ul style="list-style-type: none"> • Control Unit • Foot Control • Handpiece Contra-Angle, 20:1, Optic, X-SG20L • Micromotor with cable, 2 m, Optic, SGL70M • Sterile irrigation tubing (5pc) 	Y1500399
	Handpiece Contra-Angle, 20:1, Optic, X-SG20L	C1003-STM
	Handpiece Straight, 1:1, Optic, X-SG65L	H1009-STM
	Handpiece Straight, 1:1, Non-Optic, SGS-ES	H264-STM
	Handpiece Angled, 1:1, Non-Optic, SGA-ES	H263-STM
	Sterile irrigation tubing (5pc)	Y900113-STM
	Micromotor with cable, 2 m, Optic, SGL80M	Y1500507



International Headquarters

Institut Straumann AG
 Peter Merian-Weg 12
 CH-4002 Basel, Switzerland
 Phone +41 (0)61 965 11 11
 Fax +41 (0)61 965 11 01
www.straumann.com

Legal manufacturer

NAKANISHI INC.
 700 Shimohinata, Kanuma
 JP-Tochigi 322-8666
 Phone +81 289 64 3380
 Fax +81 289 62 5636