

PROSTHETICS

STRAUMANN®

METAL SCANBODIES

Durability and precision you can rely on.



STRAUMANN® METAL SCANBODY DESIGNED TO ENHANCE YOUR DIGITAL WORKFLOW.

Straumann's new Metal Scanbodies feature durable titanium construction and careful design considerations to improve accuracy and support a precise digital workflow. Improve your digital impressions, enhance patient comfort, and streamline your workflow.

1 ROUNDED EDGES


No blackout in CAD software.

2 ONE FLAT SIDE

Helps define position in CAD software.

3 HIGH AND CONTOURED EMERGENCE PROFILE

Reduces pressure on gingiva and optimizes visibility and scanability.

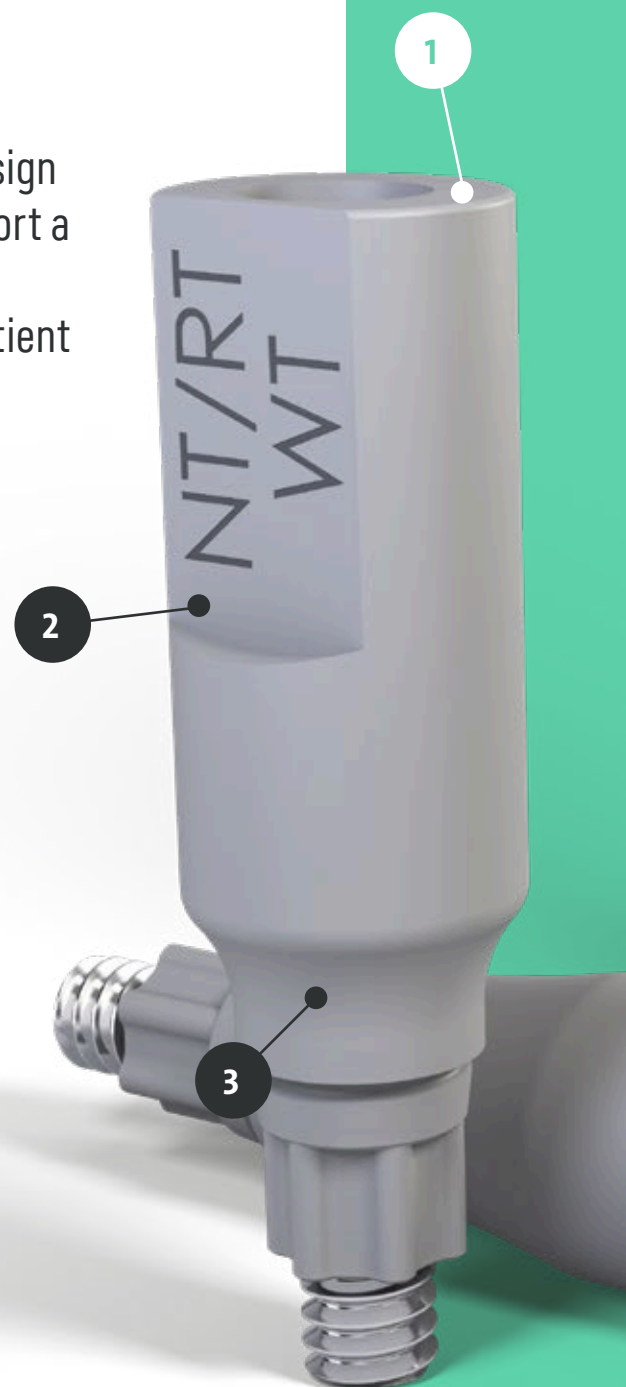



4 NARROW BODY

More space for adjacent teeth.

5 UNIQUE NON-REFLECTIVE COATING

Eliminates need for scan spray or powder.

6 LASER-MARKING

Easy to identify parts.

7 DURABLE TITANIUM


Radiopaque. Re-usable.

STRAUMANN® METAL SCANBODY

Quick reference

Implant Type	Platform	OLD PEEK Mono Scanbody	NEW Metal Scanbody
CrossFit® – Compatible to BLT Implant System	SC	025.0025	025.0077
	NC	025.2915	025.0078
	RC	025.4915	025.0079
CI RD – Compatible to Pure Ceramic Implant System	CI RD 4.1	032.041	032.130
SRA – Compatible to Screw-retained Abutments	SRA 3.5	025.0000	025.0080
	SRA 4.6	025.0001	025.0081
synOcta® – Compatible to Tissue Level Implant System	NNC	048.173	048.383
	RN	048.168	048.384
	WN	048.169	048.385
TorcFit™ BL – Compatible to BLC and BLX Implant Systems	RB/WB	065.0035	065.0103
TorcFit™ TL – Compatible to TLC and TLX Implant Systems	NT/RT/WT	036.3220	036.3230

DOWNLOAD
CAD/CAM LIBRARY

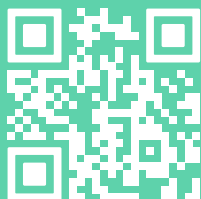

ORDER ON
E-SHOP



490.880/en/B/00_06/24

International Headquarters

Institut Straumann AG
Peter Merian-Weg 12
CH-4002 Basel, Switzerland
Phone +41 (0)61 965 11 11
Fax +41 (0)61 965 11 01
www.straumann.com

© Institut Straumann AG, 2023. All rights reserved.

Straumann® and/or other trademarks and logos from Straumann® mentioned herein are the trademarks or registered trademarks of Straumann Holding AG and/or its affiliates.

 **straumann**